
Solio Alfa⁺

- EN INSTRUCTIONS FOR USE
- DE BEDIENUNGSANLEITUNG
- FR INSTRUCTIONS D'UTILISATION
- NL GEBRUIKSIJNSTRUCTIES
- ES INSTRUCCIONES DE USO
- IT MANUALE OPERATIVO
- RU ИНСТРУКЦИИ ПО ПРИМЕНЕНИЮ

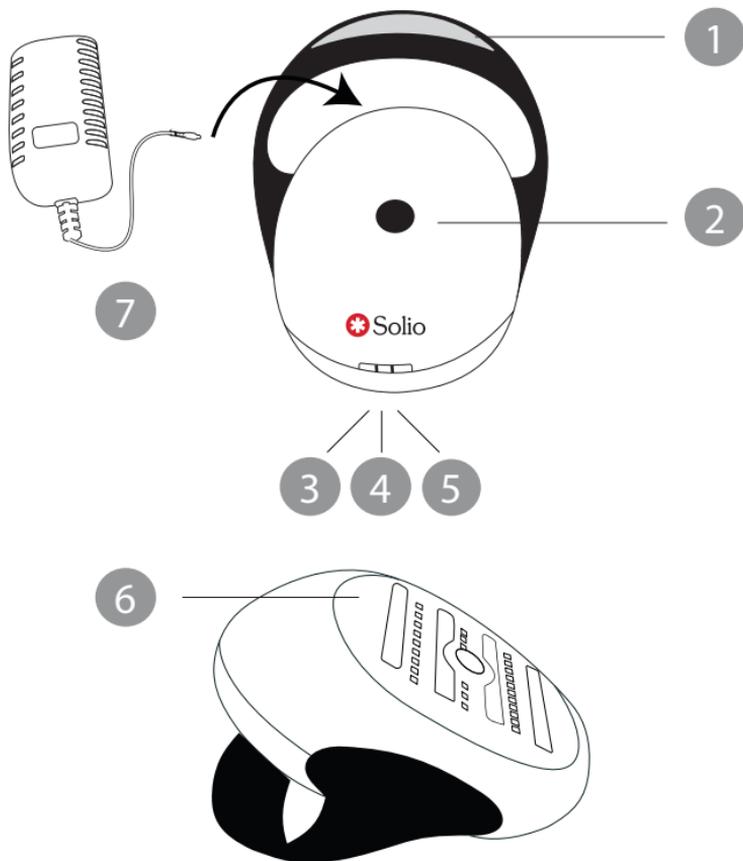Copyright 2019 © DMT Ltd.

All rights reserved.

EN

Legal Notice

DMT Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by DMT Ltd is believed to be accurate and reliable at the time of publication. However, DMT Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of DMT Ltd. Data is subject to change without notification. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from DMT Ltd. Specifications are subject to change without notice.

Solio and the Solio logo are registered trademarks of DMT Ltd, Maskit 8A, Herzliya 4673308, Israel.

Table of Contents

1. Get to know your Solio Alfa+	3
2. Contraindications & warnings	4
3. Do not use Solio Alfa+ if:.....	7
4. What to expect when using your Solio Alfa+.....	8
5. Get to know your Solio Alfa+	9
6. Get started!.....	12
7. After treatment	12
8. Maintenance and storage.....	13
9. Customer service.....	13
10. Warranty	13
11. Troubleshooting	16
12. Disposal	17
13. Specifications	18
14. Signs and symbols	19
15. Output Data.....	20
16. Labeling	20
17. Test Results	21

1. Get to know your Solio Alfa+

Solio Alfa+ is an energy-based device that provides relief from pain resulting from sports injuries, chronic pain, trauma or menstruation. It is a powerful electrical device and should be used with special attention to safety. **Be sure to read all warnings, contraindications and safety precautions before use, and follow them strictly.** We recommend you re-familiarise yourself with the instructions before each treatment.

1.1 Package contents

Your Solio Alfa+ package includes:

- Solio Alfa+ device
- Adapter with cable
- User manual
- Slider Gel
- Quick guide

1.2 Device description

Solio Alfa+ is a hand held device that provides pain relief. The device has an adjustable strap (#1), a *control button* (#2), 3 *indicator lights* (#3,4,5), a *treatment surface* with 4 *electrodes* (#6), and an *AC/DC adapter* with cable (#7). Refer to the diagram at the front of this manual to see the location of each feature of your Solio Alfa+.

1.3 Intended use

Solio Alfa+ is intended to provide pain relief. Typical uses are for the relief of aching muscles, muscle spasms, trapped nerves, sprains, stiff or swollen joints, minor arthritis, carpal tunnel syndrome, tennis elbow, tendon injuries, neuropathic pain, general inflammations, or menstrual pain.

2. Contraindications & warnings

This device must not be used by the following people: children, adults with limited physical, sensory or psychological capacities, or those lacking experience or knowledge in how to use the device. Children should be supervised to ensure that they do not play with the device.

Do not use Solio Alfa+ for any purpose other than as described in this manual.

Do not use Solio Alfa+ if you have:

- ! a history of skin cancer or areas of potential skin malignancies.
- ! applied topical steroids on the treatment area within the last month.
- ! received radiation therapy or chemotherapy treatments within the past 3 months.
- ! a pacemaker, heart arrhythmias or any other severe known heart disorder.
- ! an implantable metal device, or pierced jewelry in the treatment area
- ! a history of diseases stimulated by heat, such as recurrent herpes simplex in the treatment area, unless you have consulted your physician, and received preventative treatment before using the device.

- ! a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- ! uncontrolled diabetes.
- ! blood clotting disorders. Consult with your physician before treating.
- ! recently treated the intended area with any of the following topical substances: Alpha- Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A®, retinoids, azelaic acid, or any other physician prescribed drug. Consult your physician before use!
- ! any skin condition such as active herpes, including brown or black spots, birth marks, moles, warts, eczema, psoriasis or lesions in the treatment area.
- ! open wounds in the treatment area. Wait for the area to heal before using the device on that specific part of the body.
- ! silicone implants in the treatment area.
- ! been treated with Accutane® (isotretinoin) recently. Consult with your physician.
- ! any other condition which, in the opinion of your physician would make it unsafe for you to use Solio Alfa+.

Additional warnings

-  Do not use when you are pregnant or nursing.
-  Do not use on the face.
-  If you have any discomfort when using the default treatment mode, set the device to the treatment mode designed for small areas (2 lights on).
-  Do not use the default mode (3 lights) or small areas mode (2 lights) for the treatment of hot, swollen areas, areas diagnosed as inflamed, or areas with varicose veins. When treating these areas always set the device to inflammation mode (1 light).
-  Remove all jewelry before treatments.
-  Do not use near bathtubs, showers, basins or other vessels containing water. Do not place or store the device where it can fall, be pushed or placed anywhere near water. There is a risk of electrocution. Unplug immediately if the device comes into contact with water.
-  Never leave the device unattended when plugged into an outlet. Keep the AC/DC adapter accessible at all times, so it can be removed from the socket easily in case of emergency.
-  Do not carry the device by the AC/DC adapter cable.

-  Do not attempt to modify, open or repair your device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself.

3. Do not use Solio Alfa+ if:

- the adapter, attachments or accessories for this device are not approved by DMT Ltd.
- the device has a damaged cable or adapter. Keep the power cable away from heated surfaces.
- you see or smell smoke when it is in use.
- the device casing or treatment surface is cracked or damaged.
- the device has been dropped or has fallen onto a hard surface.
- the device is not working properly, or appears to be damaged in any other way.

In such cases, contact Customer Service.

4. What to expect when using your Solio Alfa+

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. However every energy-based therapy device, including those designed for home use, involves some degree of risk. Therefore it is important that you understand and accept the risks and complications that may occur.

During a treatment session it is normal to experience or feel:

- **a sensation of warmth:** This is typical when using the device in 'default mode' (the 3 indicator lights are on). The sensation should disappear by the end of the treatment, or within an hour. If the sensation is uncomfortable for you, continue the treatment using the setting with 2 indicator lights (refer to the section [5.2 Treatment modes](#)).
- **nothing at all:** This is typical when using the device in 'inflammation mode' (just 1 indicator light is on). Don't worry! The device still emits pain relief energy, even though you don't feel any warmth.
- **some redness:** It is not uncommon to see some mild, pink-like skin colouration during or just after a treatment. Again, this is normal and should fade within an hour or two.
- **slight swelling in the treated area:** This should also disappear within an hour or two.
- ! Stop use immediately if any of these effects are more severe.
- ! See your doctor if the effects do not go away within 2 to 3 days.

On rare occasions treating with Solio Alfa+ might cause:

- excessive swelling and redness, or burns and blisters. If this happens, immediately use a damp cool cloth to cool the area for at least 10 minutes and apply an anti-burn cream. In most cases these reactions subside within a few days.
- red spots in the treatment area which may last for a few days. If this happens, we recommend that you cool the area with chilled water immediately.
- **!** Consult your physician if these symptoms persist.

5. Get to know your Solio Alfa+

Solio Alfa+ is simple to use.

5.1 Device controls

- **Control button** (#2): Press to turn the device on, and long press to turn it off.
- **Lights** (#3, 4, 5):
All 3 lights go on when the device is turned on. Short press to cycle between the modes (refer to section [5.2 Treatment modes](#)).

Solio Alfa+ only emits energy when it is on, you have applied Slider Gel, the *treatment surface* is in good contact with the skin, and you are moving the device across your skin. **When Solio Alfa+ emits energy there will be a constant red glow from the treatment surface.**

Solio Alfa+ Slider Gel (or any other conductive gel) enables conductivity. Without the gel, the device will not work. Make sure that you always have a generous amount of gel on your skin.

5.2 Treatment modes

For maximum comfort, flexibility and effect, Solio Alfa+ has three treatment modes, to suit the type of pain and the part of the body being treated.

Types of pain and treatment area

Most pain types, and most parts of the body

Small areas (e.g, palms, ankles)

Hot, swollen areas, or areas diagnosed as inflamed, or areas with varicose veins

Lights


When the device is in 'inflammation mode' (only one indicator light is on), you will not sense the warmth that is typical when using the other treatment modes of Solio Alfa+.

5.3 Adjusting the strap

If the strap is too tight, you can remove the rubber section (#1).

- To release the section, press the top of the strap until the section pops out.
- To put the section back, line it up from below and push it back into place.

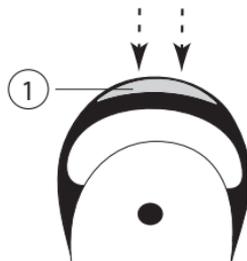5.4 Treatment Plan

A full Solio Alfa+ treatment session takes just 10-20 minutes.

Treatment intervals

- Start treating up to **twice a day** for up to **10 days**.
- If the pain persists after 10 days, continue treating **once a day** for up to **90 days**.
- If the pain persists after that, continue treating **2-3 times a week**, until the pain disappears.

Treatment Intervals

- Use Solio Alfa+ only as long as you feel pain.

6. Get started!

To use Solio Alfa+, follow these steps:

1. Remove the device and other components from box.
2. Plug the *adapter cable* into the device socket and plug the *adapter* into an electrical outlet.
3. Make sure your skin is clean and dry, then apply a thin layer of Slider Gel to the treatment area.
4. Press the *control button*. All three *indicator lights* will go on, and your device is ready.
5. Adjust the treatment mode for the area you are treating (refer to section [5.2 Treatment modes](#)).
6. Place the *treatment surface* on your skin, and move the device in a circular motion across your skin for the recommended period (refer to section [5.4 Treatment Plan](#)).

7. After treatment

When you have finished your treatment session:

1. To turn the device off, long press the *control button*.
2. Unplug the *adapter* from the electrical outlet.
3. We recommend you clean the treatment area to remove any residue of Slider Gel.

8. Maintenance and storage

It is important to clean your device after each treatment.

-  Unplug the device.
-  Use a clean damp cloth to clean the device and *treatment surface*. You may use cleaning fluids with up to 70% alcohol, but never use acetone or other abrasive liquids, because they will damage your device.
-  After cleaning, we recommend you store the device in its original box.
-  Never immerse your device or any of its parts in water.

9. Customer service

For more information, see www.soliotherapy.com. This manual is also available as a PDF download from the websites. Contact Customer Service to report unexpected operations or events, or any other problem with your device.

Email: contact@soliotherapy.com

Telephone: +972-(0)-775050300

10. Warranty

What Does This Warranty Cover?

DMT Ltd, its applicable affiliates (hereinafter: DMT Ltd) or the authorized dealer from whom you purchased your product (hereinafter: the Device) shall repair or replace the Device, if any material defects caused by DMT Ltd are discovered in the Device during the warranty period, whether in the body of the Device or its operation.

Who is entitled to warranty by virtue of this warranty certificate?

Only you, the original buyer of the Device who purchased the Device from DMT Ltd or its authorized dealers, are entitled to receive the warranty specified in this warranty certificate. In order to be entitled to the warranty service, you need to provide a proof that you are the original buyer of the Device by presenting an original receipt for purchase of the Device. This warranty is not an international warranty. The warranty will only be honored in the country in which the authorized dealer or the DMT Ltd entity from whom you purchased the Device, are located.

What period does this warranty certificate cover?

The warranty covers repair requests received within one year (1) from the Device's original purchase date. Repair requests after one year (1) from the date of purchase are not covered by this certificate of warranty. Repair requests regarding a Device that was repaired or replaced must be submitted no later than: (a) one year (1) from the original purchase date of the Device or (b) ninety (90) days from the date that the repaired or replaced Device was returned to you, the later to occur. If a sticker is attached to your Device stating a warranty period longer than that indicated above, the warranty period mentioned on the sticker will apply, and the Device will be repaired or replaced as part of the warranty until such period expires.

If you purchased the Device in an EU country, you have an additional year (according to the EU-Directive 1999/44) to claim for deficiencies or malfunctioning of your Device, which were present at the time of purchase, but were not apparent then, as determined by DMT Ltd or the authorized dealer. In order to make this extended warranty effective you must contact the authorized dealer from whom you purchased the Device or the local service center and notify them about the production fault immediately, and no later than 60 days after detecting the fault.

How to get service for the warranty?

When asking for customer services you may be asked to fill in a form for returning the product and you shall receive a reference number. In case of a defect that is not covered by the warranty for whatever reason, DMT Ltd can charge you for the shipping and repair costs if you agreed to repair the Device at your expense.

What is not covered by this certificate of warranty?

This warranty does not cover any damage resulting from an accident, misuse (including using the Device other than as described in the use instructions accompanying the Device), defective use, lack of reasonable care, irregular and unreasonable use including without limitation commercial or professional use, adding parts not supplied together with the Device, losing parts, using unauthorized spare parts or exposing the Device to any electrical current other than the indicated current. The warranty will be cancelled if the Device was treated by anyone not authorized by DMT Ltd and DMT Ltd will not be liable for any damage resulting from repairs performed by anyone not authorized by DMT Ltd. The warranty does not include replacing used non-durable parts.

WARRANTY DISCLAIMER

THE EXPRESS WARRANTY SET FORTH IN THIS DOCUMENT IS IN LIEU OF ALL OTHER WARRANTIES INCLUDING THE IMPLIED WARRANTIES OF MERCHANTABILITY, AND FITNESS FOR A PARTICULAR PURPOSE. THIS WARRANTY IS THE SOLE AND EXCLUSIVE REMEDY FOR ANY DEFECT OF THE DEVICE.

LIMITATIONS

IN NO EVENT SHALL DMT Ltd, ITS AFFILIATES OR AUTHORIZED DEALERS, BE LIABLE FOR ANY INDIRECT, SPECIAL, PUNITIVE, EXEMPLARY, INCIDENTAL OR CONSEQUENTIAL DAMAGES (INCLUDING, WITHOUT LIMITATION, DAMAGES FOR LOSS OF USE OR CHARGES FOR SUBSTITUTE EQUIPMENT), REGARDLESS OF THE LEGAL THEORY ON WHICH THE CLAIM IS BASED EVEN IF THEY HAVE BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. IN THE EVENT THAT THE REMEDIES SET FORTH IN THIS WARRANTY ARE HELD TO FAIL OF THEIR ESSENTIAL PURPOSE OR IF DMT Ltd, ITS AFFILIATES OR AUTHORIZED DEALERS ARE HELD

TO BE LIABLE FOR DAMAGES ON THE BASIS OF ANY LEGAL THEORY, IN NO EVENT SHALL DMT Ltd, ITS AFFILIATES OR AUTHORIZED DEALERS, LIABILITY FOR DAMAGES EXCEED THE WHOLESALE PRICE OF A SINGLE UNIT OF THE DEVICE AS OF THE DATE THAT THE CAUSE OF ACTION ACCRUES. YOU ASSUME ALL RISK AND LIABILITY FOR LOSS, DAMAGE OR INJURY TO YOU AND YOUR PROPERTY AND TO OTHERS AND THEIR PROPERTY ARISING OUT OF THE USE OR MISUSE OF, OR INABILITY TO USE, THE DEVICE EXCEPT TO THE EXTENT CAUSED DIRECTLY BY THE NEGLIGENCE OF DMT Ltd.

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS WHICH VARY FROM COUNTRY TO COUNTRY. SOME COUNTRIES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES OR IMPLIED WARRANTIES, SO THE ABOVE DISCLAIMERS AND LIMITATIONS MAY NOT APPLY TO YOU.

11. Troubleshooting

Do not attempt to open or repair your device. You may exposure yourself to dangerous electrical components and may cause serious harm to yourself. It will also void the warranty.

Problem	Checks
My Solio Alfa+ does not start.	Make sure the <i>adapter</i> is properly connected to the device, that it is plugged into an electrical outlet on the wall, and that you have pressed the <i>control button</i> .
There is no energy emitted when the device touches my skin.	Check that: <ul style="list-style-type: none"> • the device is connected and turned on • you have applied sufficient gel • the treatment surface is in good contact with the skin, and • you are moving the device across the treatment area

Problem	Checks
One or more lights blink continuously.	The device has malfunctioned. Unplug the device and restart.

For problems or queries, please contact Customer Service.

12. Disposal

Help protect the environment. Do not dispose with household waste. Be sure to follow your country's laws for the disposal of electrical items.


13. Specifications

Model no.	H2101/H2102
Applied part	Type BF
Protection rating	IP30 (EN60529)
Light source	Light Emitting Diodes Peak emission wavelength: 630nm \pm 10nm IR LED light wavelength: 850nm \pm 50nm
Radio Frequency	1 MHz \pm 30%; 24W max
Package size	22cm X 22cm X 14cm
System weight	200gr
Transport & storage between uses and storage condition	Temperature: minus 40° - 70°C Relative humidity: 10% - 90% Atmospheric pressure: 500 - 1060 hPa
Operating conditions	Temperature: 10° - 35°C Relative humidity: 30% - 75% Atmospheric pressure: 700 - 1060 hPa
Adapter models and rating	Model: KSAPVO361200300D5
	Input: 100-240V; 50-60Hz; 0.9A Output: 12V DC; 3A

14. Signs and symbols

	Products marked with this symbol fulfil the requirements of the EU directives.
	Degree of protection against electric shock. BF applied part.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.
	Keep dry.
	Follow the instructions for use.
	This equipment intentionally supplies non-ionizing RF energy.
IP30	Not waterproof.

15. Output Data

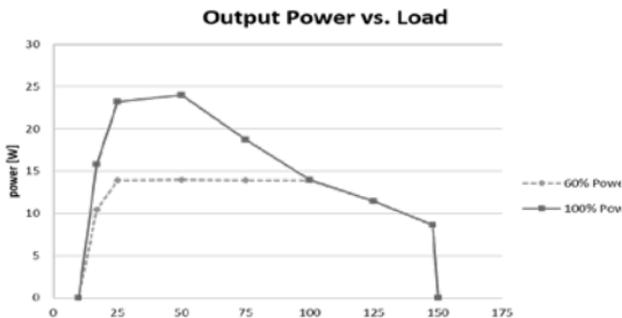16. Labeling

Protection rating against the ingress of water: Ordinary.

This device is not suitable for use in the presence of flammable anaesthetic mixtures, or with air, oxygen or nitrous oxide.

17. Test Results

Summary of Test Results

Test	Standard	Class/ Severity Level	Test Result
Documentation (IEC 60601-1-2 sections 4 and 5)			
General requirements General requirements for EMC.	section 4.1	--	Complies
External labels	section 5.1	--	Complies
Conformity of Users' Manual	section 5.2.1	--	Complies
Accuracy of Technical Description	section 5.2.2	--	Complies
Emission (IEC 60601-1-2 section 7)			
Conducted emission Freq. range:150 kHz - 30 MHz	CISPR 11	Group 1 Class B 230, 120 & 100 VAC mains	Complies
Radiated emission Freq. range: 30 - 1000 MHz	CISPR 11	Group 1 Class B	Complies
Harmonic current emission test	IEC 61000-3-2	230 VAC mains	Test NA
Voltage changes, Voltage fluctuations and Flicker test	IEC 61000-3-3	230 VAC mains	Complies
Immunity (IEC 60601-1-2 section 8)			
Immunity from Electrostatic discharge (ESD)	IEC 61000-4-2	8 kV contact discharges & 15 kV air discharges	Complies
Immunity from radiated electromagnetic fields	IEC 61000-4-3	10.0 V/m; 80 MHz + 2.7 GHz, 80% AM, 1 kHz	Complies
Immunity from Proximity field from wireless communications equipment	IEC 61000-4-3	List of frequencies from 9 V/m up to 28 V/m, PM (18 Hz or 217 Hz), FM 1 kHz	Complies
Immunity from Electrical Fast transient (EFT)	IEC 61000-4-4	± 2.0 kV; Tr/Tf - 5/50 ns, 100 kHz on AC mains	Complies
Immunity from Surge	IEC 61000-4-5	±1.0 kV DM; Tr/Tf - 1.2/50 (8/20) µs on AC mains	Complies
Immunity from conducted disturbances induced by radio- frequency fields	IEC 61000-4-6	3.0 VRMS & 6.0 VRMS on AC mains; 0.15+ 80 MHz, 80% AM 1 kHz	Complies
Immunity from voltage dips, short interruptions and voltage variations	IEC 61000-4-11	0 % - 0.5 cycle & 1 cycle; 70% - 25 cycles; 0% - 250 cycles on 230 & 100 AC mains	Complies

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